




RIVAL™

○ REDUCE

□ VIEW

▶ BITE

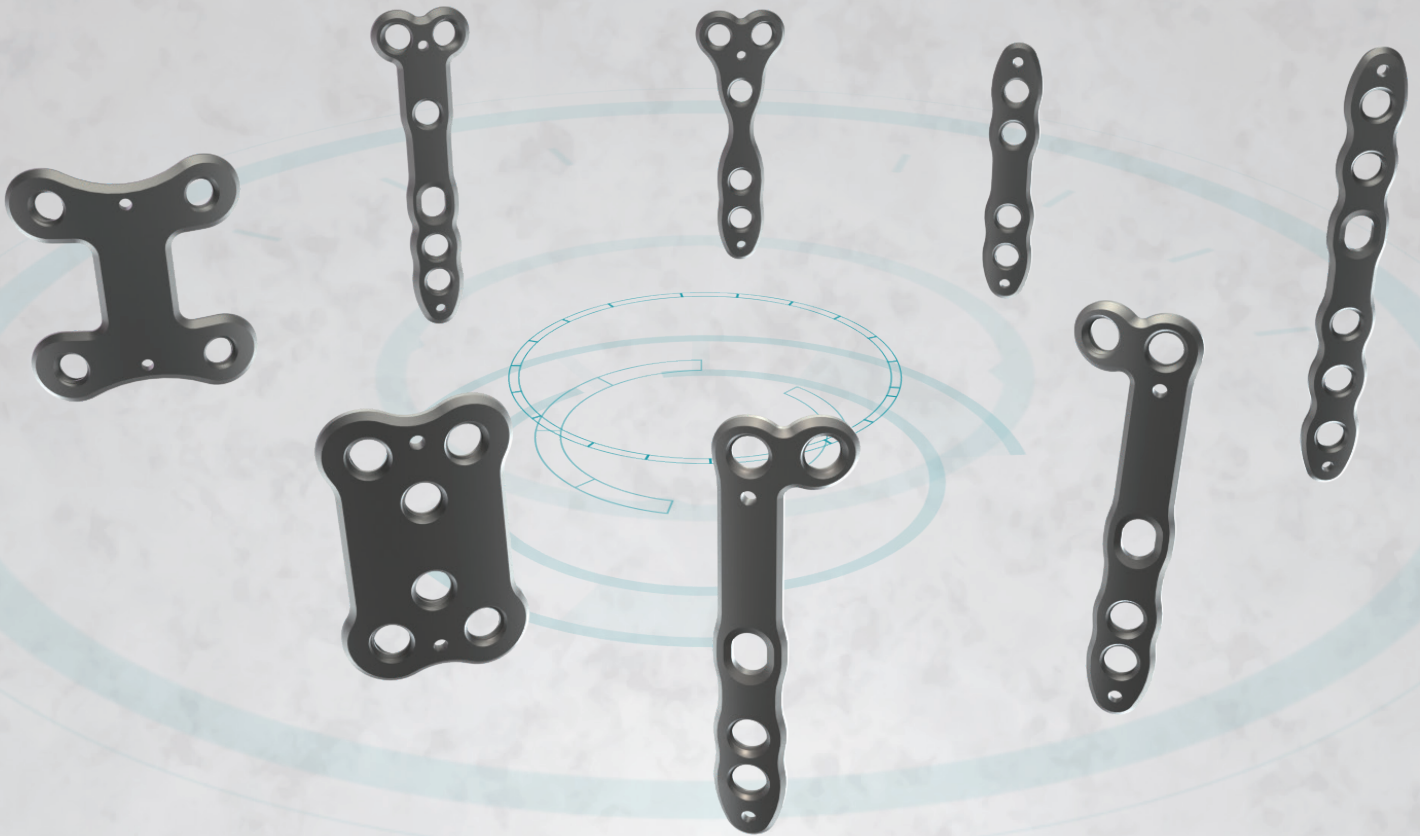

○ REDUCE

Fracture Plating System

Wide range of implants and sizes


Optimal implant selection

Designed to

reduce soft tissue irritation


Designed to

reduce screw prominence

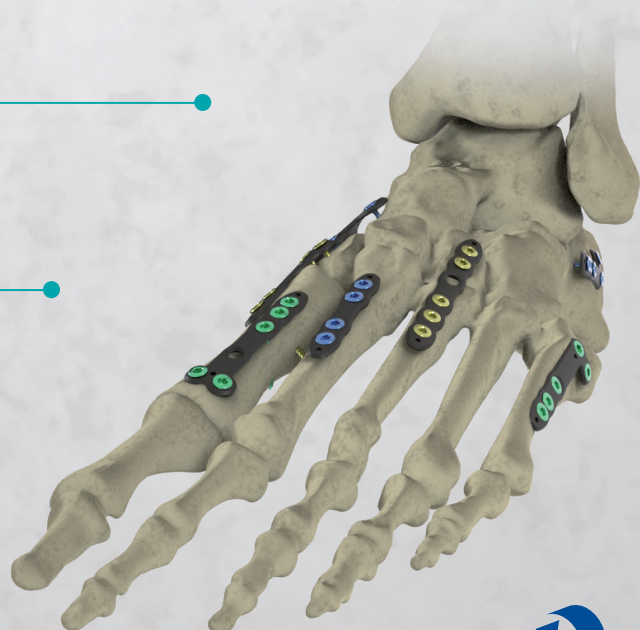


Designed to be used with 2.4mm, 2.7mm and 3.2mm locking and non-locking screws

Locking screws

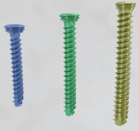


Non-locking screws

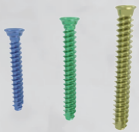

Double lead thread screw design

For fast and controlled insertion

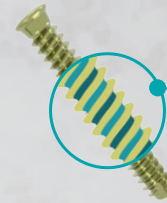

Oblong compression hole

Intraoperative compression up to 1.1mm with the corresponding non-locking screws

Guide wire holes

For easy implant placement during surgery

Type II anodised titanium for **Increased fatigue resistance**



Low profile

1.3mm plate thickness



Screws lock flush into plates reducing screw prominence




Universal instrumentation

Designed


- to be used across all RIVAL™ product lines
- to help with inventory management
- to provide a comprehensive range of instruments

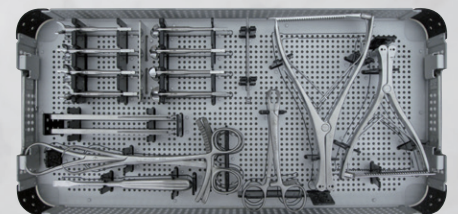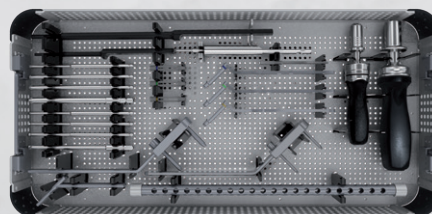

Distributed by:

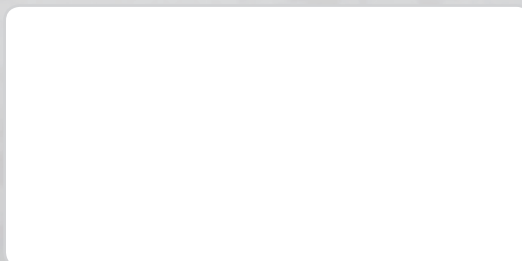

Manufactured by:
 ORTHOFIX Srl
 Via Delle Nazioni 9, 37012 Bussolengo (Verona), Italy
 Telephone +39 045 6719000, Fax +39 045 6719380

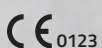

www.orthofix.com

RIRE-1904-PL-E0 A 10/19

