


TL-HEX COMPONENTS

Full Rings

Part Number	Description
56-20000	Full Ring, 140mm, TL-HEX
56-20020	Full Ring, 160mm, TL-HEX
56-20040	Full Ring, 180mm, TL-HEX
56-20060	Full Ring, 200mm, TL-HEX
99-56-20080	Full Ring, 220mm, TL-HEX (Sterile)

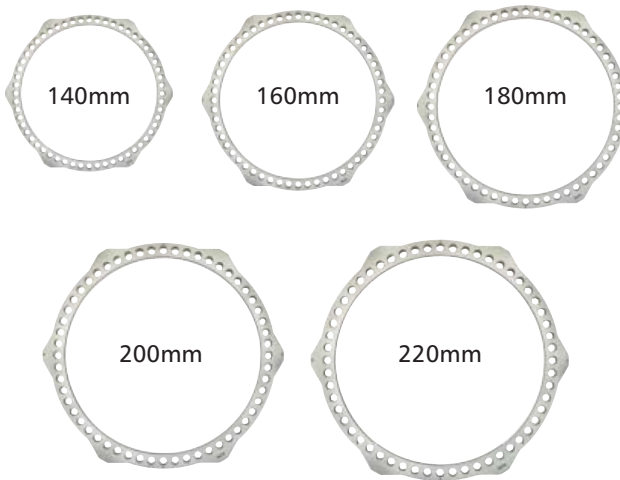

Struts

Part Number	Length
50-10200	Short Strut TL-HEX - 92mm-122mm
50-10300	Medium Strut TL-HEX - 114mm-184mm
50-10400	Long Strut TL-HEX - 158mm-318mm

Short



Medium



Long




TL-HEX Strut Exchange Kit
cod. 50-10500

Part Number	Description	Quantity
50-10180	TL Plus Medium Quick Adjust Strut	1
50-10190	TL Plus Long Quick Adjust Strut	1
52-10210	TL+ AO 1/8" Hex Power Driver Attachment	1
52-1018	TL - 1/8" Hex Driver	1
54-1154	TL - Wrench, Combo, 10mm	1
54-2226	TL - 90 Degree Tubular Wrench	1
50-13020	TL+ Extended Nut M6 X 1	2
50-10214	Strut Direction Clip Kit, TL-HEX (6 pcs)	1

5/8 Rings

Part Number	Description
56-21000	5/8 Ring, 140mm, TL-HEX
56-21020	5/8 Ring, 160mm, TL-HEX
56-21040	5/8 Ring, 180mm, TL-HEX
56-21060	5/8 Ring, 200mm, TL-HEX
99-56-21080	5/8 Ring, 220mm, TL-HEX (Sterile)

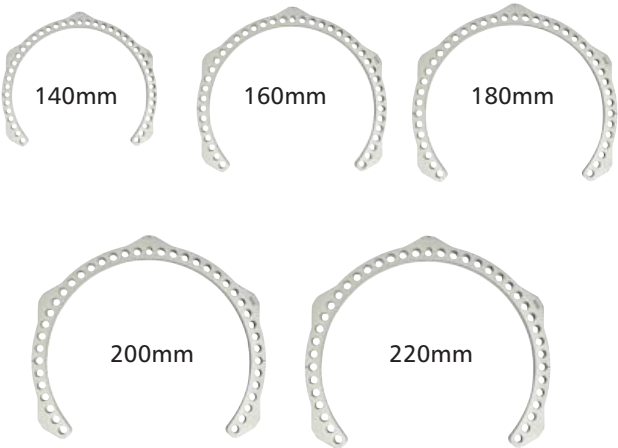

Number Clips

Part Number	Description
50-10215	Strut Number Clip Kit, TL-HEX (6 pcs) Consisting of: Strut Number Clip, #1, TL-HEX Strut Number Clip, #2, TL-HEX Strut Number Clip, #3, TL-HEX Strut Number Clip, #4, TL-HEX Strut Number Clip, #5, TL-HEX Strut Number Clip, #6, TL-HEX



Direction Clips

Part Number	Description
50-10214	Strut Direction Clip Kit, TL-HEX (6 pcs)






TL-HEX TRAYS

Rings Tray

The Ring Tray can accommodate a combination of the TL-HEX rings for sterilization. The maximum number of rings that the tray can hold is four full rings and two 5/8 rings.

Part Number	Description
56-20000	Full Ring, 140mm, TL-HEX
56-20020	Full Ring, 160mm, TL-HEX
56-20040	Full Ring, 180mm, TL-HEX
56-20060	Full Ring, 200mm, TL-HEX
56-21000	5/8 Ring, 140mm, TL-HEX
56-21020	5/8 Ring, 160mm, TL-HEX
56-21040	5/8 Ring, 180mm, TL-HEX
56-21060	5/8 Ring, 200mm, TL-HEX
30110086	Rings Tray, TL-HEX (empty)

Material out of the tray	
99-56-20080	Full Ring, 220mm, TL-HEX (Sterile)
99-56-21080	5/8 Ring, 220mm, TL-HEX (Sterile)

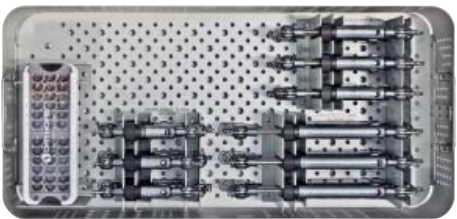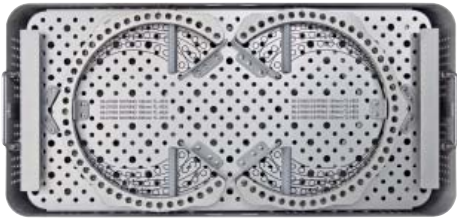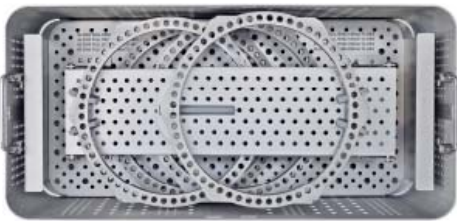

Struts Tray

Part Number	Description	Quantity
50-10214	Strut Direction Clip Kit, TL-HEX	3 Kits (18 Strut Direction Clips)
50-10215	Strut Number Clip Kit, TL-HEX	3 Kits (3 of each Strut Number Clip)
50-10200	Short Strut, TL-HEX - 92mm-122mm	6
50-10300	Medium Strut, TL-HEX - 114mm-184mm	6
50-10400	Long Strut, TL-HEX - 158mm-318mm	6
30110085	Struts Tray, TL-HEX (empty)	

HEX...cellent
COMPATIBILITY

All TL-HEX components are fully compatible with the TrueLok Ring Fixation System. Please refer to: TL-1201-OPT and TL-1204-PL.

Manufactured by:
ORTHOFIX Srl
Via Delle Nazioni 9
37012 Bussolengo (Verona)
Italy


Telephone +39 045 6719000
Fax +39 045 6719380



Your Distributor is:



www.tlhex.com

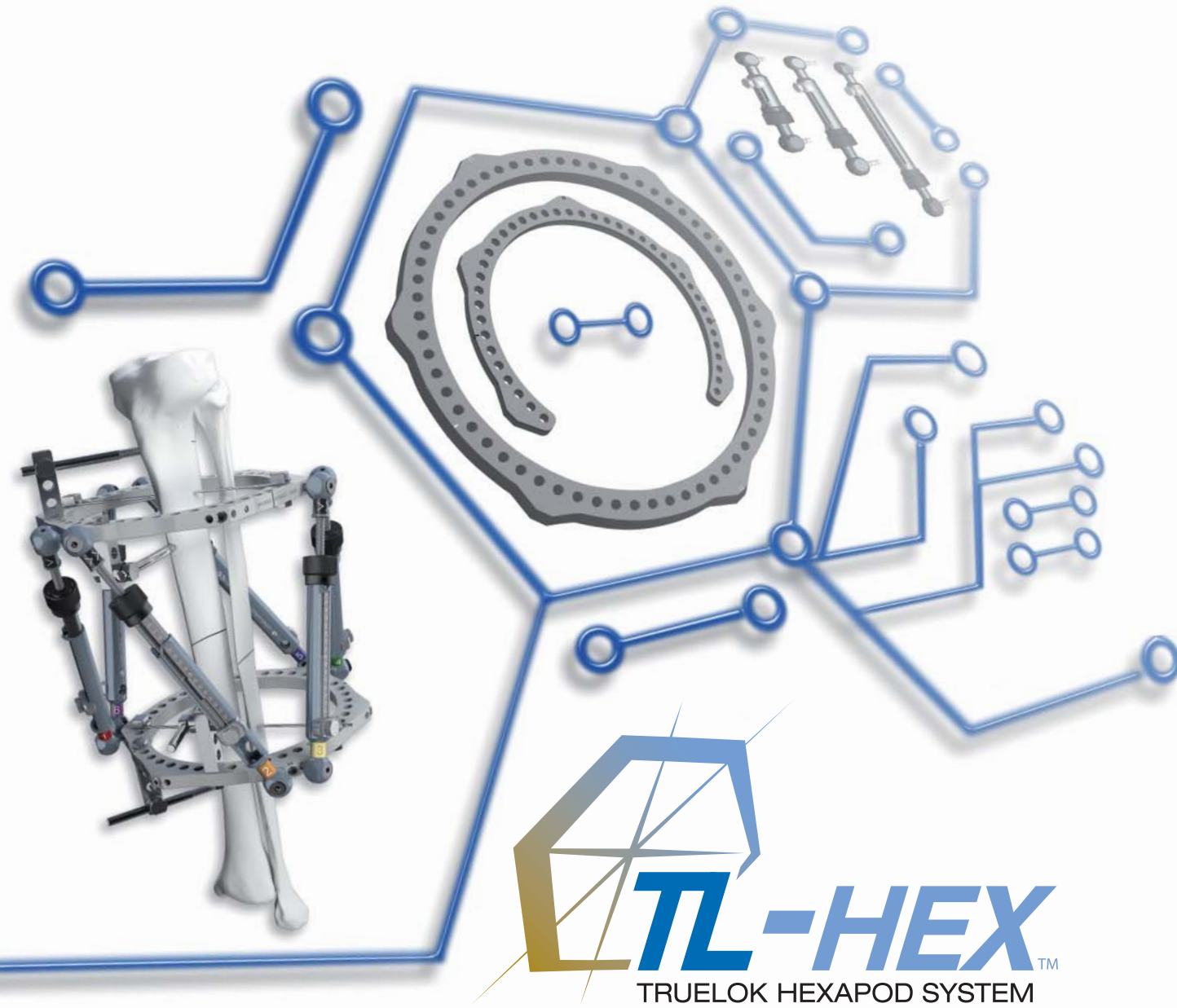

The future of the hexapod
system built on solid circular
fixation experience

TL-HEX...cellent!

HEX...cellent EXPERIENCE

TL-HEX and Truelok Ring Fixation System were designed in collaboration with Texas Scottish Rite Hospital for Children and preserve most of the original principles and methodology of Professor Ilizarov, while providing modern hardware and software

TL-HEX provides exceptional stability due to its unique

Stable

aluminum-stainless steel and metal-plastic interface

TL-HEX provides simplified Hardware and

Simple

Software for both Deformity and Trauma management

The distinctive Strut design allows performing acute and gradual

Versatile

adjustment in deformity correction and complex trauma procedures

HEX...cellent PERFORMANCE

EXCELLENT STABILITY
Unique strut head design increases frame stability

FAST LOCKING
One set screw locks down both struts



RAPID ADJUSTMENT
TL-HEX Struts consist of two telescoping aluminum tubes, which can be locked together at various lengths using the side locking bolt and clamp washer

EASY ADJUSTMENT
Simple pull and click method for patient adjustment

HARDWARE



STRUTS

TL-HEX struts are available in three sizes (short, medium, long), providing adjustment range from 92mm to a maximum extension of 318mm

RINGS


TL-HEX full rings and 5/8 rings are available in 5 sizes with internal diameters of 140mm, 160mm, 180mm, 200mm and 220mm



HEX...cellent ADJUSTMENT

CLIPS

The **number clips** indicate the strut number and the **direction clips** indicate the direction of gradual adjustment to simplify the patient interaction with the frame

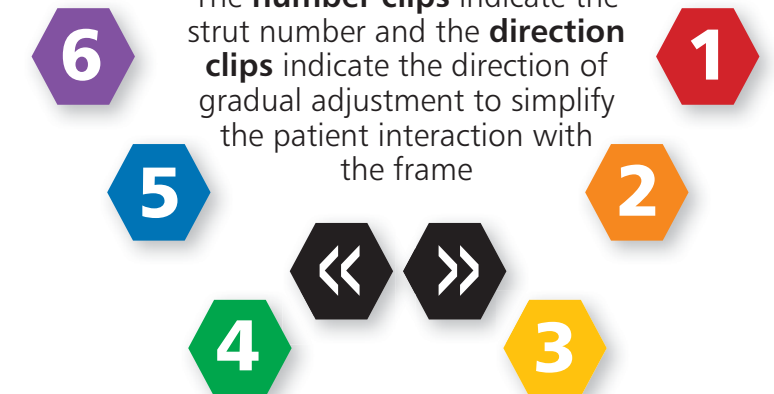

HEX...cellent CONTROL

SOFTWARE

The intuitive and user friendly web-based TL-HEX software is empowering all surgeons who want a comprehensive support throughout pre/intra/post operative phases

